

## FIXATION SET, ORTHOPEDIC [TOC](#)

STRYKER/HOWMEDICA

MODEL# 4920-9-952

Manufacturer's PHONE: (201) 831-5726 FAX: (201) 831-6553

Website:

<http://www.trauma.stryker.com/products/hoffmann2.php?tab=&vis=&loc=>


AAC "L" ITEM NSN 6515-01-546-0559

UA's: 267B, 623A, F302, F308, N302, N308, N342, N542, N583, N742, N783

USAMMA POC: Customer Relations Management, 301-619-4301/4316, DSN 343-4301/4316

UNIT PRICE: \$2,002.00

**System Description:** External fixation is a method of immobilizing bones to allow a fracture to heal. External fixation is accomplished by placing pins or screws into the bone on both sides of the fracture. The pins are then secured together outside the skin with clamps and rods. The clamps and rods are known as the "external frame." Some of the advantages of external fixation are that it is quickly and easily applied. The risk of infection at the site of the fracture is minimal, but there is a risk of infection where the pins are inserted from the skin into the bone. The Forward Deployed Field Pack is for temporary stabilization and transportation only.

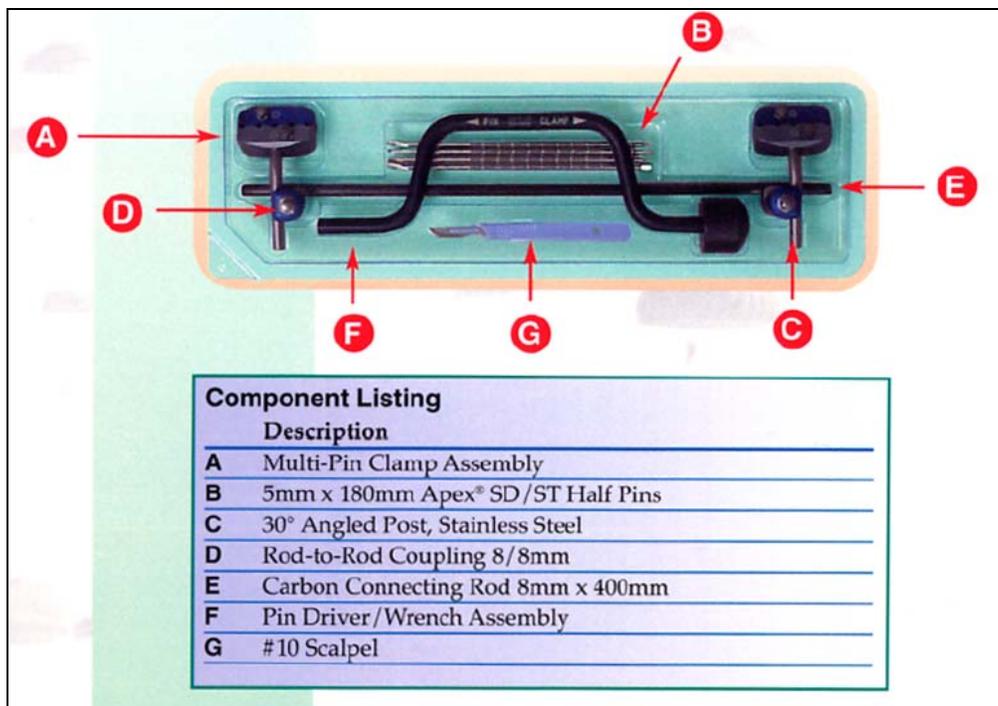

FIXATION SET, ORTHOPEDIC

MODEL# 4920-9-952

NSN 6515-01-546-0559

**FIXATION SET, ORTHOPEDIC  
PHYSICAL & ENVIRONMENTAL SPECIFICATIONS**

<b>Dimensions</b>	<b>Height</b>	<b>16.5 inches</b>
	<b>Width</b>	<b>7 inches</b>
	<b>Depth</b>	<b>2 inches</b>
<b>Cube</b>		<b>0.13 cu ft</b>
<b>Weight</b>		<b>1.9 lb</b>

**FIXATION SET, ORTHOPEDIC  
(START-UP)**

NSN	NOTES	Nomenclature	PN	QTY	UI	Unit Price	Total Price	DIST/MFG	Shelf Life, months	Refrig Y/N
6515015460559		Fixation external ortho ster Hoffman II	4920-9-952	1	EA	\$2,002.00	\$2,002.00	Howmedica Osteonics	24	N
6515014631475		Pin clamp assembly	4920-2-020	2	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
6515015131261		Pin bone 5mm dia x 180 mm long	5018-6-180	7	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
6515014631497		Post angled fixation ortho 30° SS	4920-2-140	2	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
6515014631498		Coupling fixation rod-to-rod ortho 8 / 8mm	4920-1-010	2	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
6515015040434		Rod carbon connecting 400mm long x 8mm wide	5029-8-840	1	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
No NSN required		Pin driver / wrench assembly	5194-6-000	1	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N
No NSN required		# 10 scalpel		1	EA	Included /W 0559 NSN	Included /W 0559 NSN	Howmedica Osteonics	N/A	N